

尺寸: (80X115)mm

材质: 135g硫酸纸

Prospera Magic TENS PL009A OTC TENS Electronic Pulse Massager

MAGIC TENS OPERATION STEPS

1 Install



Connect two buttons on the the electrode pad to the main unit;

2 Get ready



Remove the TENS unit from its attached paper, and place it on your neck, shoulder, arm or leg area;

3 Start TENS



Press Power button (indicate the shape or log of the button) for 3 seconds, once beep sound heard, green dot light indicates the unit is ready to work;

4 Intensity Adjustment




Press **+ / -** Power button to adjust **intensity**, two "beeps" indicate at the lowest / **highest** intensity;

5 Mode Adjustment



Single-click on/off button, it will beep once each time to change to a different mode; total 5 modes;

6 Finish up

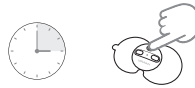

Auto Shut-off: the Magic TENS will shut off automatically in 15mins;
Manually Shut-off: press Power button (show the image) for 3 seconds to shut off.

Prospera Corporation
sales@prosperecorp.com

448 Commerce Way, Livermore, CA 94551
Tel: (925) 225 0888