



Shiatsu Foot Massager
+heat +air

DL-002

User Manual

The DL002 user manual provides detailed product specifications, operation and warning information. Please read this manual and understand it thoroughly before installation and operation. Keep this manual for future reference.

CONTENTS

Safety Precautions.....	1
Product Description.....	2
How to Operate.....	3
Troubleshooting.....	4
Care & Maintenance.....	5
Product Specifications.....	6

About Prospera Company

Prospera is a leading US provider of high quality, innovative, affordable health consumer electronics. Prospera designs and manufactures a variety of products from FDA approved TENS units (510K), massagers, foot spa to other health and beauty products. Our products help reduces muscle tension and release stress and back pain caused by today's active lifestyle. Learn more about Prospera at www.prosperacorp.com

About the DL002 Shiatsu Foot Massager +heat +air

The Prospera Shiatsu Foot Massager is an ergonomic light weight device to work on your feet exclusively for comfort and convenience.

The foot massager comes with deep tissue Shiatsu kneading massage, air compression massage and heating to relive fatigue. It helps you to relax after a long walk or a day of work.

The foot covers can be removed for washing to keep the sleeves clean and prevent smell.

Auto 15 min shut off timer built in the device.

You will feel renewed. Can you feel it?

SAFETY PRECAUTIONS


1. People with the following conditions should consult their physician before using the product:
 - Who uses a medical electronic device, which is embedded in the body, such as a pacemaker.
 - Who suffers from heart disease.
 - Who suffers from backbone injury or disease.
 - Who is suffering from bone fragility, especially on the backbone.
 - After surgery.
 - Who is pregnant or in a period of monthly illness.
 - Who suffers from malignant tumors.
 - Who is suffering from osteoporosis.
 - Who is suffering from blood circulation problems or thrombosis.
 - Who is suffering from skin disease or infection.
 - Who is undergoing medical treatment or feels some physical abnormality previously or is being treated now.
 - Who has some abnormality or deformity of limbs or backbone.
 - Who has injuries, diseases due to illness or accident.
2. Do not let water enter the product or be splashed near the surface of the product.
3. Water contact with electricity is dangerous. Do not leave the product in a humid environment (such as in the bathroom or shower room or swimming pool).
4. Before connecting the product to the power supply, check that the product voltage and supply voltage are the same.
5. If the plug, power cord or product itself is damaged, do not use the product.
6. In order to prevent the risk of electric shock, do not open the product outer cover.
7. Product is not intended for use by the physically or mentally impaired (including children), unless there is supervision or guidance by the people responsible for their safety when they use the product.
8. If the product or the power cord is damaged, do not try to repair the product. This must be handled by Prospera authorized service technicians otherwise the product guarantee will become void.
9. Children should not use this product. They should be under supervision to ensure that they will not treat this product as a toy.
10. If an unusual sound or noise is being heard from the product, turn off the product and immediately disconnect the power connection. Please contact an Prospera authorized service technicians.
11. Do not use the product if cover is damaged or if there is a mechanical or electronic failure. Contact the Prospera authorized service technicians
12. Only use the product in accordance with the instructions in this user manual.
13. Do not use this product in high temperature conditions or beside hot objects such as heaters, stoves.
14. Do not apply medical or cosmetic creams on the product or on your foot before using this product.
15. Do not destroy, bend, stretch, twist or tie the power cord. Damaged power cord can cause fire or electric shock.
16. People who are sensitive to heat should not use the heat function of product.
17. Recommended usage time for an individual is 15 minutes at a stretch.

PRODUCT DESCRIPTION

1. Control Panel

2. Foot Slot

3. Cloth Cover

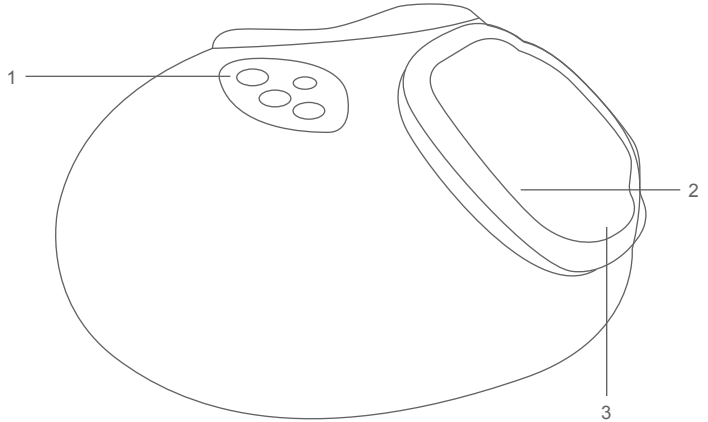

Control Panel

A. MODE

"AUTO" Mode : A Combination of Kneading and Air Massage for Feet


"MASSAGE" Mode: Kneading Massage for Soles

"AIR" Mode : Full Air Massage for Feet

B. INTENSITY levels

C. HEAT Therapy

D. STOP

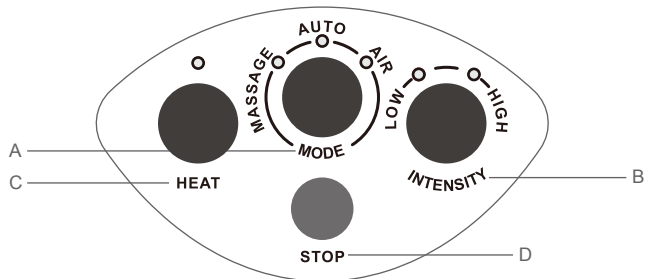

How to operate the DL002

1. Press auto mode button to turn device on default program, which has a combined kneading and air press massage.
2. Press auto one more time to shift to deep tissue kneading massage mode.
3. Press auto 3rd time to switch to air press massage. Use intensity button to choose low or high intensity.
4. Press heat button for heating feature.
5. Once you finish, press stop to complete the massage session.

The key functions are as follows:

No.	Button	Function
1.	MODE	Press to turn ON with default "AUTO" mode of low intensity Press again to turn on "MASSAGE" mode Press once more to turn on "AIR" mode Press again to stop all massages
2.	STOP	Press to stop the massage any time
3.	INTENSITY	2 Levels of intensity: HIGH / LOW (Applicable to "AUTO" and "AIR" modes only)
4.	HEAT	Heat therapy of the soles

Auto Timer: 15 minutes. Power will turn Off automatically after 15 minutes to prevent over heating.

Always unplug the AC power cord after use.

TROUBLE-SHOOTING

This page summarizes the most common problems you could encounter with the product. If you are unable to solve the problem with the information below, please contact prosperacorp.com

Problem	Possible Cause	Suggestion
The product does not work	The product is not plugged in properly.	Ensure the voltage on the product corresponds to the local mains voltage and then plug in again.
	The product has just finished its auto timer session and has automatically switched to standby mode.	If you want another session, press the MODE button again.
	You may have used the product continuously for 30 minutes or more where the over heat protection has been activated.	Unplug the product and let it cool down for 60 minutes before you use it again.
	An excessive load has been applied on the product.	Reduce the excessive load and start again.
The product is abnormally warm	The product is likely to be overused.	Unplug the product and let it cool down before using it again.
Noise can be heard from the product	Noise is made by the mechanical parts of the product.	This is normal and does not mean that the product is faulty.

CARE & MAINTENANCE

Cleaning

Do not use steel wool, abrasive cleaning agents or corrosive liquids to clean the product.

1. Unplug the plug from the power outlet.
2. Wipe with a dry cloth to clean the product.
3. Do not use a damp cloth or detergent.

Storage

1. Unplug the plug from the power outlet.
2. Store the product in a cool dry place away from sunlight.
3. Do not put a heavy weight load on the product.

If the cover is damaged or there is any mechanical failure, please stop using the product.

PRODUCT SPECIFICATIONS

Name of Product	Prospera Shiatsu Foot Massager
Model Number	DL002
Heating	Yes
Power Input	AC 100-240V 50-60Hz/0.8A
Power Consumption	18W
Auto Timer	15 mins
Product Dimensions	13"x13"x6.5"
Box Dimensions	14.3"x13.8"x6.9"
Nett / Gross Weight	Approx. 3.15lbs/4.0lbs
Safety Compliance	CE and UL

Consumer Limited Manufacturer's Warranty

Within the period of this warranty, Prospera will repair or replace, free of charge, any part proving defective in material and workmanship. All warranty repairs and service must be performed by Prospera authorized technician and service facility.

Warranty exclusions:

This warranty does not apply to any costs, repairs, or services for the following:

1. Repairs necessitated by use other than normal home use.
2. Damaged resulting from misuse, abuse, accidents, alternation or improper installation.
3. Corrective work necessitated by repairs made by anyone other than Prospera authorized service technicians.

Non-warranty service is provided on a "per incident" basis. Consumer will verify that the product has failed and provide information for servicing the unit. All applicable repairs, parts, shipping, handling, local tax and "per incident" fees will be charged for non-warranty repairs and support calls.

For warranty service, please call 925-292-5978, M-F, 8:00am to 4:00pm, Pacific Standard Time. Or send email to sales@prosperacorp.com. Consumer may also mail service request to this address:

Prospera Corporation
Warranty Service
448 Commerce Way, Suite A
Livermore, CA 94551

Warranty card

Warranty Registration Form

Product Name: _____

Product Model: _____ Purchase Price: _____

Date of Purchase: _____ Retailer: _____

First Name: _____ Last Name: _____

Address: _____

City: _____ State: _____ Zip: _____

Phone: _____ Email: _____

Signature and date: _____

Please mail to: sales@prosperacorp.com

Tel: (925) 292 5978

Prospera Corporation
448 Commerce Way
Livermore, CA 94551
www.prosperacorp.com